



Right on the nose

Respironics Simplicity Nasal Mask bridges the gap between performance and comfort.

Sometimes achieving effective therapy can compromise comfort. So we developed Simplicity to eliminate pressure points, while still delivering high performance.

Simplicity is a small, lightweight nasal mask worn just on the tip of the nose. By having no contact with the bridge of the nose, it helps avoid common pressure points. As a result, both men and women wearers can enjoy a positive seal without irritation. Patients also will enjoy an unobstructed view – in fact, they can even wear their eyeglasses. Simplicity is one more way Philips Respironics is providing sensible solutions that help make it easier for you to fit 100% of your patients, 100% of the time.

PHILIPS

RESPIRONICS


The exhalation port results in quiet air flow, and its innovative built-in design eliminates the need for additional valves or ports.

360° swivels at the top and bottom allow tubing to twist freely as the patient moves.

The Deluxe Softcap headgear is a lightweight, streamlined mesh cap with adjustable straps that allow for optimal fit.



The Simplicity headgear is specially designed for a secure, comfortable fit. Made of soft, flexible material, it features a unique five-strap design that evenly distributes tension, adds stability, and helps to prevent mask leaks. Velcro tabs allow easy adjustment for a custom fit.

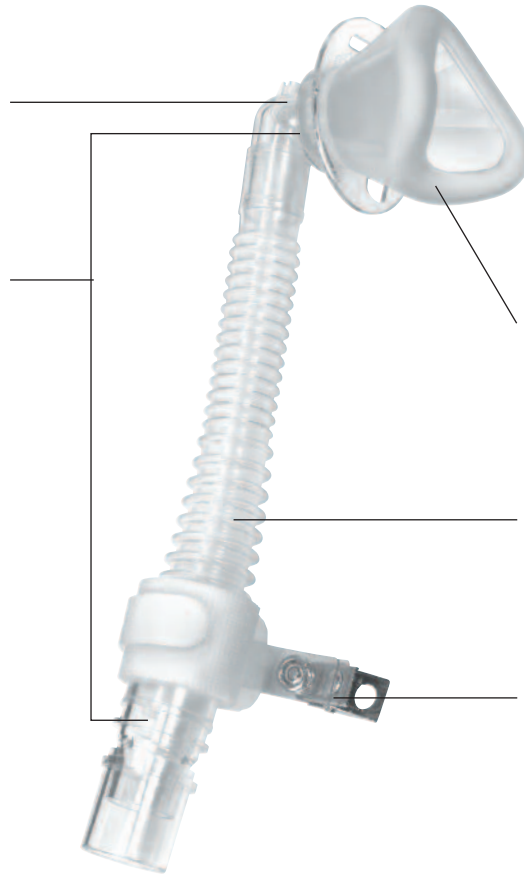

The mask cushion is made of soft, pliable silicone that gently contours around the tip of the nose to provide a snug seal.

Flexible, lightweight tubing is permanently attached to the mask and connects easily to standard circuit tubing.

The Stability Clip on the tubing may be attached to sleepwear or the headgear for added support. A Velcro tab secures placement of the clip on the tubing.

Item	Part number
Simplicity mask	
Mask with Simplicity headgear, small	1002761
Mask with Simplicity headgear, medium	1002763
Parts and accessories	Part number
Simplicity headgear	1002765
Small replacement cushion	1003681
Medium replacement cushion	1003680
Tubing assembly with Stability Clip	1004209

Respironics, Simplicity, and Softcap are trademarks of Respironics, Inc. and its affiliates, all rights reserved. Velcro is a trademark of Velcro Industries.



© 2009 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Broudy KB 10/19/09 MCI 4103059 PN 1025657

Philips Healthcare is part of Royal Philips Electronics

Respironics Asia Pacific
+65 6298 1088
Respironics Australia
+61 (2) 9666 4444
Respironics Europe, Middle East, Africa
+33 1 47 52 30 00
Respironics United Kingdom
+44 800 1300 845

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service
+1 724 387 4000
800 345 6443 (toll free, US only)
www.philips.com/respironics