



Unsurpassed Quality

DeVilbiss offers the longest standard warranty in the industry on our IntelliPAP CPAP Systems, up to five years of coverage, thanks to our two-year extended warranty program (3 + 2).

- Rest easy knowing your products are covered for the full life of care
- Less than 1% return rate over five-year period
- 100% tested to meet or exceed our strict quality and sound level standards



Cost-free Compliance

Get the data you need without the cost. Through on-board SmartCode® Technology, healthcare providers and clinicians can obtain critical therapy information quickly.¹


- No add-on charges
- Enter the SmartCode data into the report generator at www.IntelliPAP.com
- A report will be immediately generated that includes usage and clinical efficacy² information




Usage Data				
	Last 90 Days	Last 30 Days	Last 7 Days	Last Day
% Days of at least 4 hours (not continuous) in a 24 hour period	64%	90%	100%	100%
Day Count	90	30	7	1
Days at least 4 hours	58	27	7	1
AHI	8.5	5.5	5.0	4.5
High Leak Flow Time	5%	2%	3%	1%
NRI	1.75	1.50	1.00	0.50
EPI	2	2	3	1
SmartCode	7MTJT-4594-7DMAQ	D0Y4-3D97-OPNE	3CCC-379D-7CNF	F34F-383-3PRD
While Breathing Hours	5.4			
Adherence Score (4 continuous hours)	67% (20/30)			
Titration Mode	CPAP	CPAP	CPAP	CPAP
Minimum Usage Threshold	4 Hours	4 Hours	4 Hours	4 Hours

1. Available on IntelliPAP® Standard, Standard Plus, AutoAdjust®, Bilevel, and AutoBilevel Systems.
2. IntelliPAP Standard Plus, AutoAdjust, and AutoBilevel will report therapy efficacy data.

Complete Comfort

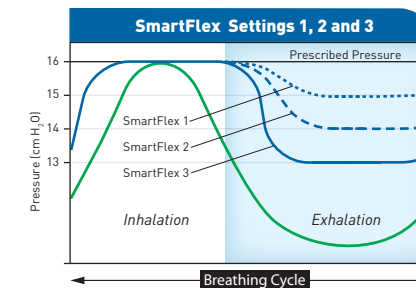

SmartFlex™ Technology
SmartFlex eases pressure during exhalation to allow patients to exhale against lower pressures.³

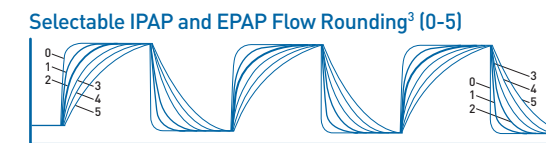

Flow Rounding™ Technology
Patented Flow Rounding provides an adjustable easing of delivered flow for inhalations and exhalations.⁴

3. U.S. Patent 5,865,173; Application No. 5545595 filed on 11/06/1995.

Humidification

The Heated Humidification System provides optimal humidification and additional comfort for your sleep therapy. The integrated heated humidifier seamlessly works with each IntelliPAP model.

Compact Design



With one of the smallest footprints in the industry, the IntelliPAP's unique stacked design – even with humidifier – takes up minimal space on a nightstand.

- Whisper-quiet operation
- No bulky power brick
- Ultimate travel CPAP



3. Available on IntelliPAP Standard Plus and AutoAdjust Systems.
4. Available on IntelliPAP Bilevel S and AutoBilevel Systems.